



REBOUND® POST-OP KNEE

Easy to fit. Easy to use.




REBOUND® POST-OP KNEE. EASY TO FIT. EASY TO USE.

Rebound Post-Op Knee is designed with two priorities in mind: ease of fitting for practitioners and ease of use for patients.

Simple and intuitive.

Smart-Fit™ packaging and a tool-free design make fitting fast and easy, while a numbered frame and color-coded patient touch points simplify patient education and product use. Just tell your patients, “Blue is You™.”

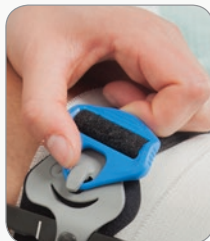

Total control.

Use the OneTouch™ hinge to quickly and easily set patient range-of-motion. Then, secure your range-of-motion settings in place with lockout clips to further reduce the risk of patient non-compliance.

Comfortable and secure.

Up to 30% lighter than the competition, Rebound Post-op Knee has ergonomic paddles and contouring struts that hug the anatomy to promote user comfort and an Anti-Migration System (AMS) wrap to help ensure secure suspension.

30%
LIGHTER*

Something for everyone.

With universal and pediatric versions suitable for most users, plus left/right contoured models for when a more anatomical fit is required, you have everything you need to fit any patient you encounter.



Contoured



Universal



Pediatric



PRODUCT HIGHLIGHTS

- 1 Large buckles featuring hook closures and grip-friendly ridges for simple brace application and removal
- 2 Numbered frame for patient ease-of-use
- 3 Color-coded drop lock and buckles for easy patient education and product use – Blue is You™
- 4 Lock-out ROM pull tabs offer -10° to 120° of ROM control, with OneTouch™ drop-lock from -10° to 30°
- 5 Anti-Migration System wrap for secure suspension
- 6 Quick-fit buckle with a lever arm helps ensure optimal strap tension and brace suspension
- 7 Telescoping frame features easy-to-engage OneTouch™ button with secure locking feature
- 8 Alligator strap tabs enable easy strap trimming
- 9 Moldable, anatomically contoured struts help ensure a precise fit



SMART-FIT™ PACKAGING

Rebound Post-op Knee features Smart-Fit™ packaging to aid in proper brace positioning and fit. Simply open the box, lay it flat, slide it underneath the patient's leg and follow the fitting instructions printed on the box.



CONTINUUM OF CARE

Össur offers a breadth of Injury Solutions that follow a recommended continuum of care model. After knee injury or surgery, protect and treat by restricting range-of-motion with the Rebound® Post-Op Knee, while managing pain with Össur cold therapy solutions. In some cases, rehab with a Functional Healing solution may be necessary. Finally, return to pre-injury activity with CTi® or Rebound® DUAL. Choose Össur products as rehabilitation progresses to follow the Össur Continuum of Care.

1 POST-OP



Rebound® Post-Op Knee

Controlled range-of-motion in a light-weight, quick-fitting design.



Cold Rush®

Convenient, whisper quiet and easy-to-use, Cold Rush delivers hours of cold therapy when and where you need it.

2 REHAB



Rebound® PCL

Designed to optimize rehab by reducing load on the PCL by applying a dynamic, physiologically-correct force during non-surgical or post-surgical treatment.



Rebound® Cartilage

A protective, functional solution designed to support regeneration of the knee cartilage after injury or surgery.

3 RECOVERY



CTi®

With anatomically-correct Accutrac® hinges, the rock-solid carbon frame provides a rigid exoskeleton that stabilizes the knee joint.



Rebound® DUAL

Functional support for ligament instabilities, including ligament instability with osteoarthritis, in a single, versatile aluminum brace.

Visit ossur.com/rebound for more information.

FOLLOW ÖSSUR ON



USA (800) 233-6263
WWW.OSSUR.COM

CANADA (800) 663-5982
WWW.OSSUR.CA

