Koalaa

How to fit the CLM

A step-by-step guide to fitting the CLM for users with above-elbow limb differences or amputations.

Let’s get started!

How to fit the CLM

This guide is here to help you confidently fit a Koalaa sleeve. Inside, you’ll find clear, step-by-step instructions on how to correctly fit the CLM to a user, attach and swap out tools, and make simple adjustments to ensure a comfortable fit.

Need help? Reach out to our support team

At Koalaa, we’re committed to making it easy for you to get the support you need. Whether you have questions about measurements, fitting, ordering, or anything else related to our products and services, our team is here for you. Just get in touch, we’re always here to help.

hi@yourkoalaa.com

+44 (0) 7482 165 098

@KoalaaCommunity

www.yourkoalaa.com

> **These tools open up opportunities to do things that matter, things that bring happiness!**

Scan to access the Koalaa Clinic Portal

For information on fitting a CLM including guides and instructional videos, visit our Clinic Portal to find everything you need.

What you need to fit a CLM

Before you get started, make sure you have everything you need to fit the CLM and follow the steps in this guide for fitting each component. The table below gives you a quick overview of each component and what it's used for.

CLM Blank

Provides the main support around the upper arm, includes the harness and inner foam to create a secure and comfortable fit

CLM Dock (Elbow)

Connects the CLM Blank to the CLM Elbow Joint + Forearm and tools

CLM Elbow Joint + Forearm

Provides movement at the elbow and supports the forearm section of the CLM

CLM Dock (Wrist)

Connects the Koalaa Tools to the CLM Forearm, allowing you to easily attach and change tools

CLM Cover

Covers the upper arm and forearm sections, providing additional support while giving the CLM a finished look

Koalaa Tools

A range of functional tools that can be attached to the CLM to support different activities and everyday tasks

What you will find in your CLM Fitting Kitt

4x CLM Blanks
Small, Medium, Large, X-Large

5x CLM Docks
2x Small, Medium, Large, X-Large

4x CLM Covers
Small, Medium, Large, X-Large

6x CLM Tools
Janet, Kitty, Amy, Rushton, Sarah, Nicole


6x Tool Foams

Fitting Equipment


Measurements required for fitting the CLM

CLM Blank Guide

Use the measurements below to determine the size of the CLM Blank

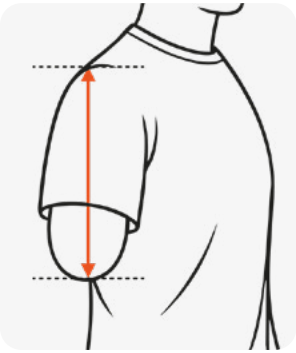Measurement G
Length from the acromion to end of limb

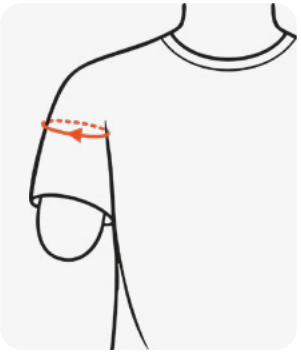Measurement I
Circumference of the top of limb

G	I	CLM Blank	Elbow Joint + Forearm	CLM Cover
16-29cm	14-19cm	S	S	S
16-29cm	19-25cm	M	M	M
22-37cm	25-33cm	L	L	L
22-37cm	33-41cm	XL	XL	XL

Scan for our online sizing calculator

Enter the user's limb measurements, and the calculator will determine the correct size CLM sleeve

Measurements required for fitting the CLM


CLM Elbow Dock Guide

Once you have chosen the correct CLM Blank size, use the measurements below to determine the size of the CLM Elbow Dock

Measurement J
Circumference of the end of limb

	S Elbow Dock	M Elbow Dock	L Elbow Dock	XL Elbow Dock
S CLM Blank	<15cm	<18cm	x	x
M CLM Blank	<15cm	<18cm	<23cm	x
L CLM Blank	x	x	<23cm	>23cm
XL CLM Blank	x	x	<23cm	>23cm

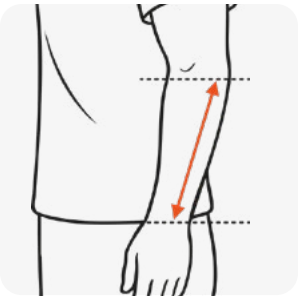CLM Wrist Dock Guide

Once you have chosen the correct CLM Elbow Joint + Forearm size, use the measurements below to determine the size of the CLM Wrist Dock

Measurement S
Length from elbow to wrist (sound arm)

	S Wrist Dock	M Wrist Dock	L Wrsit Dock	XL Wrist Dock
S Forearm	11-31cm	X	x	x
M Forearm	25-37cm	17-25cm	11-17cm	x
L Forearm	x	29-37cm	11-29cm	x
XL Forearm	x	x	17-37cm	x

Fitting the CLM Blank


Select the CLM Blank that best matches the user's limb circumference

Place the CLM Blank onto the user's limb and check the fit

Fold back the outer cover of the CLM Blank

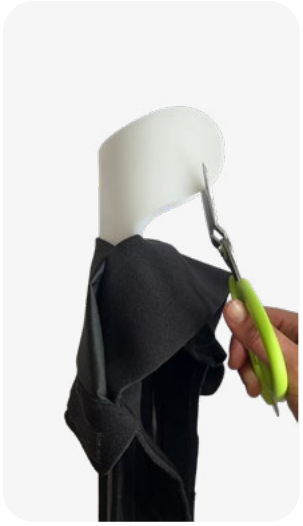Trim the wings using scissors as required to achieve an appropriate fit

If required, use a heat gun to mould the wings to desired shape

Fold the outer cover back into place

Fitting the CLM Dock (Elbow)


Select the CLM Dock to be used as the elbow unit and add the alignment guide

Position the CLM Dock on the inside of the CLM Blank

Align the CLM Dock with the user's opposite elbow

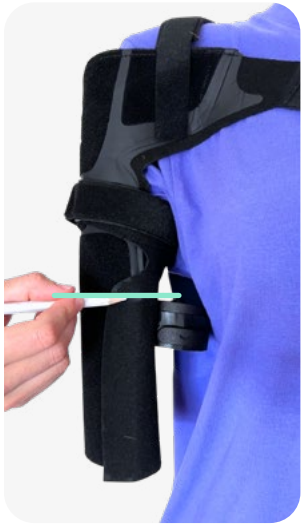Using a chalk pencil, mark the trim position at the top of the dock

Add 15 mm (1.5 cm) below this mark. This additional length will sit inside the recess of the dock

Remove the CLM Blank and trim along the marked position using scissors

Fitting the CLM Dock (Elbow)

Add the CLM Dock with the button facing forwards from the chest

Tighten all screws evenly around the CLM Dock

Ensure all screws sit flush with the surface of the dock

Add the foam inserts into the CLM Blank to achieve the correct fit for the user's limb length

Start with the orange foam to create a firm base

Add the gel foam closest to the user's limb for a comfortable, cushioned fit

Fitting the CLM Elbow Joint + Forearm

Select the correct CLM Elbow Joint + Forearm size

Attach the CLM Elbow Joint + Forearm into the CLM Dock by screwing it into position


Slide the CLM Dock inside the end of the CLM Elbow Joint + Forearm

Align with the opposite wrist and mark at the top of the CLM Dock

Add 15 mm (1.5 cm) below this mark

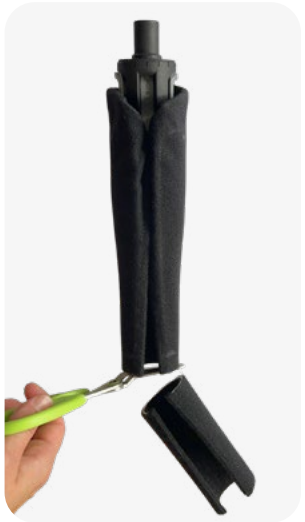Trim off the excess forearm material

Fitting the CLM Dock (Wrist)

Select the CLM Dock to be used as the wrist unit

Push it onto the end of the trimmed forearm, with button facing upwards

Tighten all screws around the CLM Dock and ensure all screws sit flush

Take the forearm section of the CLM Cover and align it with CLM Forearm, mark desired trim line


Trim off the excess cover material


Place the trimmed cover back onto the forearm

Fitting the CLM Harness

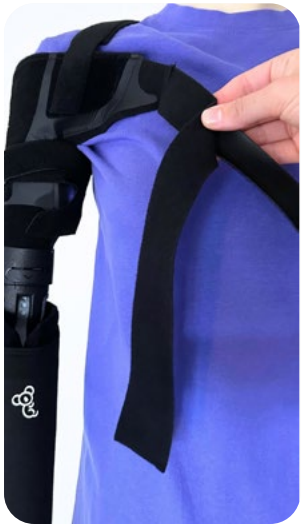Take the front strap of the harness and separate the Velcro

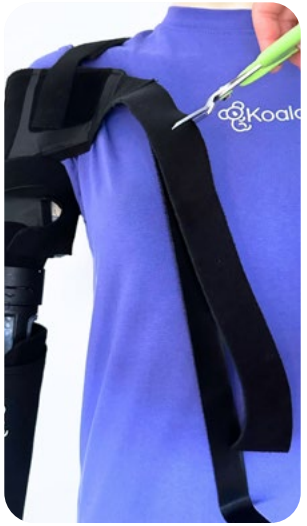
Trim the strap to the desired length so that it sits comfortably across the front of the chest

Attach the buckle strap to the trimmed harness

Take the back strap and tighten it to provide the desired level of support

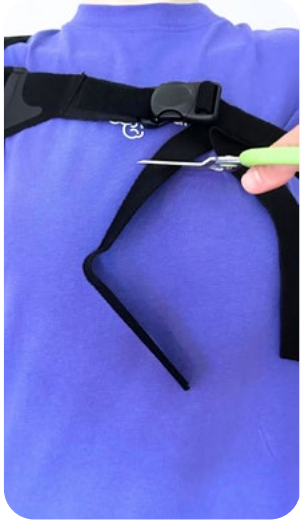
Trim off any excess Velcro and strap material

Attach the buckle and check that it is securely fastened

Fitting the CLM Harness

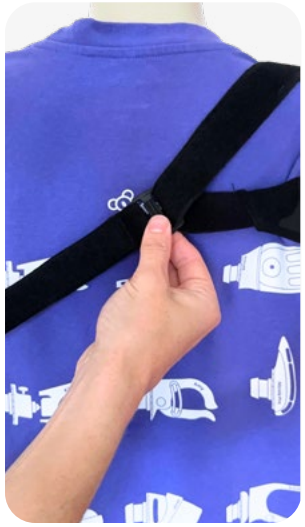Ensure the back buckle is centred on the user's back

Adjust the shoulder strap by trimming it to the required length

Secure the Velcro tab back into position on the shoulder

Adjust the upper-arm straps by loosening or tightening them as required

Trim any excess Velcro

Ensure all straps and Velcro are secure and sit flush

Fitting the CLM Cover

Take the upper section of the CLM Cover

Drape the CLM Cover over the upper arm

Fold back any excess cover material to reveal the dock and mark the required trim line


Trim off the excess cover material


Drape the cover back over the upper arm and secure the cover using the Velcro tabs

Using the CLM Elbow Joint

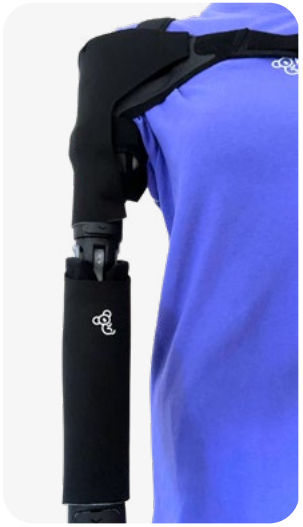You can rotate the forearm in different directions

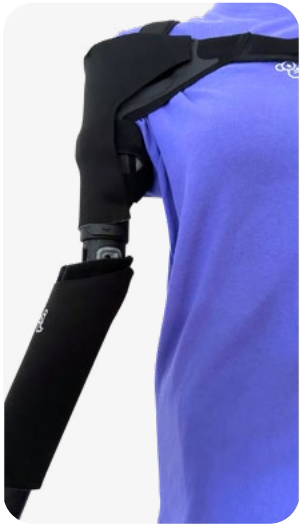Turn the forearm left or right

Lock it in position by pushing the button on the dock down

You can bend the forearm into different positions

Move the forearm up or down

Lock it in position by pressing the elbow button in

Adding tools to the CLM Dock

Select a Koalaa tool, push the tool firmly into the dock until it clicks

Rotate the tool to the desired position

Lock the tool into place by pushing the button on the dock down, ensure the lines on the tool and dock are lined up

Push the button in to release the tool

How to put on the CLM

Take the fully assembled CLM and place the user's limb into the CLM

Position the harness over the user's shoulder

Loop the chest strap under the armpit and connect the buckle at the front

Connect tools into the end of the CLM Dock at the wrist

Remove the CLM Elbow Joint + Forearm by roating it until it releases from the dock

Connect tools into the end of the CLM Dock at the elbow

CLM Max Weight Loading

Tool at Elbow: Pull / Carry

Straight Forearm: Pull / Carry

Straight Forearm: Push

Bent Forearm: Carry

Range of compatible tools

Janet

Kitty

Amy

Rushton

Sarah

Sarah Pro

Nicole

Nicole Pro

Grace

EDGA

Cosmetic

Joanie

How to care for the CLM

To keep the CLM in great condition and working at its best, a little care goes a long way. Below are some top tips to help keep it clean, comfortable.

Ensure the CLM is **comfortable** and not too tight. If the sleeve leaves marks on the skin that don't fade within 20 minutes, it's too tight and should be loosened.

Make sure the sleeve is positioned correctly, with the tool release button facing **upwards**.

Make sure all **digits** are clear of the dock when attaching tools to the CLM.

Regularly **clean** the CLM and tools by wiping away any excess dirt with a damp cloth to maintain hygiene and functionality.

To **wash** the CLM, remove the tool and inner foam, then place the CLM in a regular wash at **30°C (86°F)**. Be sure to put it in a pillowcase or laundry bag to prevent damage to other items in the machine.

To extend the life of the CLM and tools, store them in a **cool, dry place** when not in use.

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