How to use your CLM

A step-by-step guide on using the CLM for users with above-elbow limb differences or amputations.

Let’s get started!

How to use your CLM

This guide is here to help you get started with your Koalaa sleeve, so you can feel confident using it daily or for the activities you enjoy. Inside, you’ll find simple, step-by-step instructions on how to put on your CLM, how to use and swap out the tools, and how to make small adjustments for your comfort.

> **These tools open up opportunities to do things that matter, things that bring happiness!**

Need help? Your limb buddy is here for you!

At Koalaa, there’s no need for formal appointments or long waits, our limb buddies are always here and ready to help. Whether you need advice, support, or just a friendly chat, you can reach out anytime.

We’re here to help with things like tool swaps, returns, or if something just doesn’t feel quite right, just get in touch!

hi@yourkoalaa.com

+44 (0) 7482 165 098

@KoalaaCommunity

www.yourkoalaa.com

We’re with you every step of the way!

From the moment you receive your CLM, our team is here to support you on your journey. Whether you’re getting to know your new tool, learning how to use it in daily life, or adapting to something new, we’re always here to help. As you grow in confidence and begin exploring more activities, Koalaa will continue to support you, offering guidance, new options, and helping you discover our full range of tools and sleeves.

Scan Me!

Scan the QR code to access helpful how-to videos and other useful links!

How to put on your CLM

Take your CLM and place your limb into the upper sleeve

Position the shoulder harness over your shoulder

Loop the chest strap under your armpit and connect the buckle at the front

Connect the tool to the end of the forearm

You can remove the forearm by rotating it out of the dock


Connect tools to the end of the CLM sleeve

How to use the CLM Elbow

You can rotate the forearm in different directions

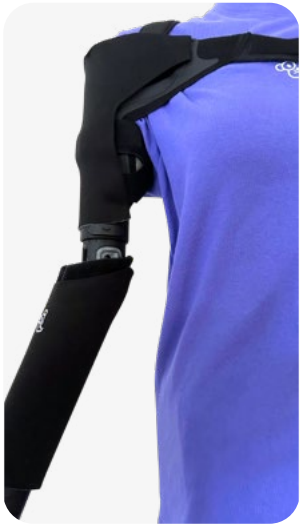Turn the forearm left or right

Lock it in position by pushing the button on the dock down

You can bend the forearm into different positions

Bend the forearm up or down

Lock it in position by pressing the elbow button in

How to add your tools

Select a Koalaa tool, push the tool firmly into the dock until it clicks

Rotate the tool to the desired position

Lock the tool into place by pushing the button on the dock down, ensure the lines on the tool and dock are lined up

Push the button in to release the tool

CLM Max Weight Loading

Tool at Elbow: Pull / Carry

Straight Forearm: Pull / Carry

Straight Forearm: Push

Bent Forearm: Carry

How to use the

Amy Tool

Max
0.5 kg
1 lbs

How to use the

Rushton Tool

Max
0.5 kg
1 lbs

How to use the


Janet Tool


Max
30 kg
66 lbs

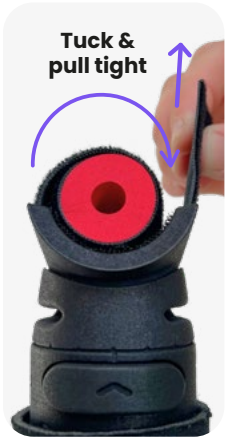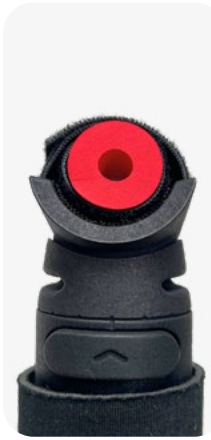How to use the


Kitty Tool


Max
30 kg
66 lbs

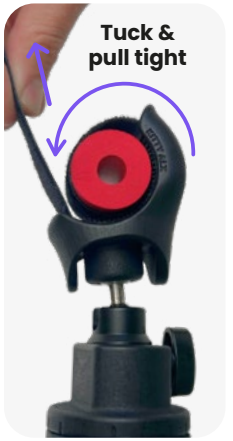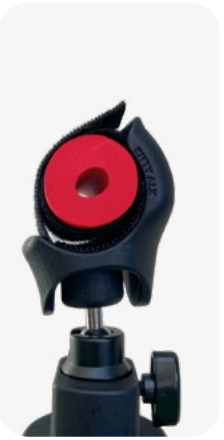How to use the

Sarah Tool

How to use the

Sarah[Pro] Tool

Max
30 kg
66 lbs

How to use the

Nicole Tool

Max
30 kg
66 lbs

How to use the

Nicole Pro Tool

Max
30 kg
66 lbs


How to use the

Grace Tool

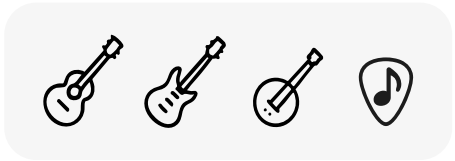
How to use the

EDGA Tool (SB&ST)

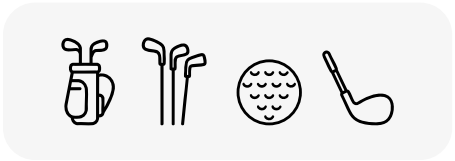How to use the functional

Cosmetic Hand

Max
0.5 kg
1 lbs

How to use the

Joanie Tool

How to care for your CLM

To keep your CLM in great condition and working at its best, a little care goes a long way. Below are some top tips to help you keep it clean, comfortable, and ready for whatever you've got planned.

Ensure the CLM is **comfortable** and not too tight on your arm. If the CLM leaves marks on your skin that don't fade within 20 minutes, it's too tight and should be loosened.

Make sure the CLM is positioned correctly, with the tool release button facing **upwards.**

Make sure all digits are **clear** of the dock when attaching tools to the CLM.

Regularly **clean** your CLM and tools by wiping away any excess dirt with a damp cloth.

To **wash** your CLM, remove the tool and inner foam, then place it in your regular wash at **30°C (86°F)**. Be sure to put it in a pillowcase or laundry bag to prevent damage to other items in the machine.

To extend the life of the CLM and tools, store them in a **cool, dry place** when not in use.

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