Cabo Elbow





INSTRUCTIONS FOR USE


Adjustable Elbow Brace


DESCRIPTION OF SYMBOLS


MD Medical Device

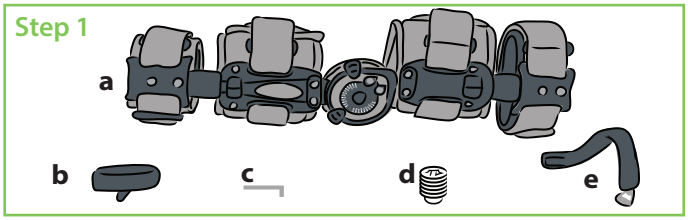


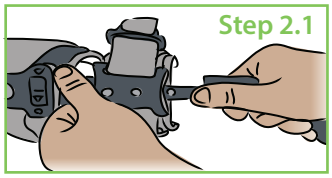


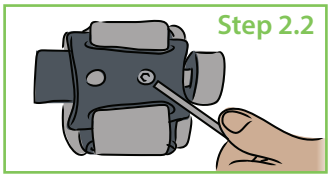


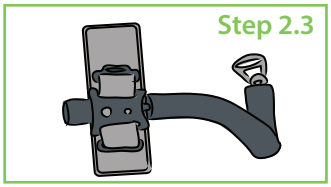


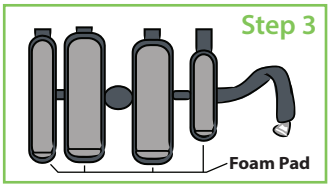


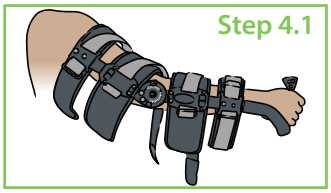


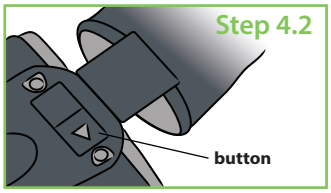


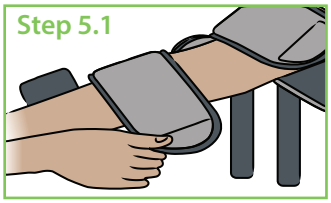


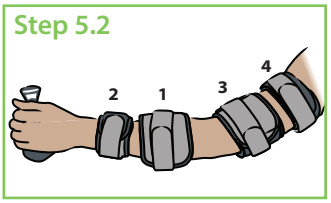


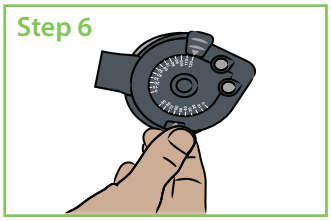


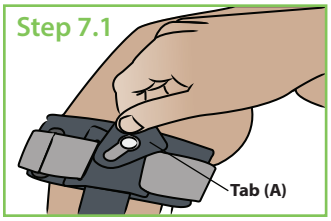

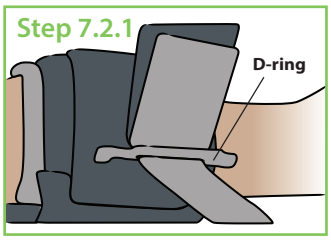

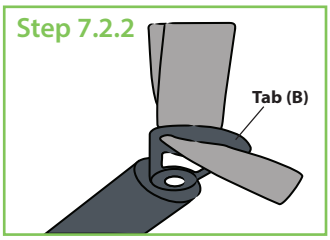

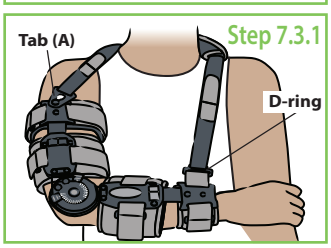

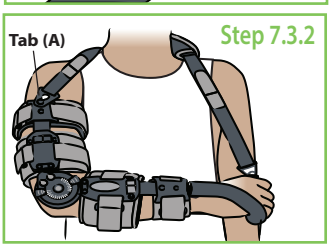

PRODUCT DESCRIPTION

The Cabo Elbow Brace is an elbow orthosis that combines an adjustable range of motion hinge and telescoping arms to insure maximum control of the elbow joint. Available in left or right configurations. Optional hand attachment for additional control.

INDICATIONS FOR USE

Intended use is to limit range of motion following surgery or injury to the elbow. Ideal for stable fractures, Biceps repair, triceps repair and strains and sprains of ligaments and soft tissue of the elbow.

GENERAL SAFETY INSTRUCTIONS / WARNINGS AND PRECAUTIONS

Read the instructions carefully before using.
Do not remove elbow brace unless instructed by your medical treatment professional. Do not change range of motion of hinge settings without permission by your medical professional. If you experienced increased pain, swelling, skin irritation, or any adverse reactions while using this product, consult your medical professional immediately.
Caution: Single patient use only product should not touch the wound directly.

ASSEMBLY/FITTING INSTRUCTIONS

Step 1: Check components before installation. Make sure all the components are inside the package: a) Brace body; b) Arm sling; c) Hex wrench; d) Hex screw; e) Hand grip assembly.

Step 2: Install the hand grip. (Hand grip optional. Consult physician.) Insert the hand grip to the hole of the forearm end cap. (Step 2.1). Install the hex screw, and use the hex wrench to tighten screw. (Step 2.2). The completed installation of the hand grip as shown in (Step 2.3).

Step 3: Open the arm cuffs.

Step 4: Hold the brace up to the arm with the hinge centered over the elbow joint. The forearm component should extend to the wrist. (Step 4.1). To adjust the length of the brace, press the button and slide the bars until they match the length of the forearm. (Step 4.2).

Step 5: **Support the injured arm:** Hold the brace next to the arm with the hinge aligned with the elbow joint. start with pad #1 on forearm closest to the hinge. Pull to desired tension and secure pad #2 the pad closest to the hand. Pull to desired tension and secure Pad #3 the pad on the upper arm closest to the elbow. pull to desired tension and secure Pad #4 the pad at the top of the upper arm.

Step 6: **Control range of motion:** Range of motion can be limited by adjusting the range of motion indicator tabs. Pull the tab up and slide to the desired setting. Let go of the tab to secure setting. Extension settings between 10° hyperextension to 120°, and flexion setting between 10° to 120°.

Step 7: **Apply shoulder strap:** Position the shoulder strap over the patient's shoulder and head with the shoulder pad located behind the neck. Connect Tab A to the end of the upper arm upright (step 7.1). If using the hand grip connect tab B to the end of the hand grip (step 7.2.2). Feed the strap through the D ring of the tab unless you are not using the hand grip. If not using the hand grip feed strap through D ring at the wrist cuff (7.2.1) . Pull to desired tension and secure hook and loop. The completed installation as shown in step 7.3.1 (without hand grip) or step 7.3.2 (with hand grip).

CLEANING AND CARE

Remove pad and hand wash with mild detergent, rinse thoroughly and let air dry. Do not use heat.

LIABILITY

ManaMed does not assume liability for the following:
- Device not maintained as instructed by the Instructions For Use.
- Device assembled with components from other manufacturers.
- Device used outside of recommended use condition, application, or environment.

SYMBOLS GLOSSARY

For an explanation of symbols, please refer to www.manamed.com/glossary/symbol-glossary.pdf

MANAMED
ManaMed, LLC
2612 Sirius Dr.
Denton, TX 76208
Ph: 888-508-0712
www.ManaMed.com
Patent No. 10,123,899

©2024 ManaMed LLC. All rights reserved.
ManaMed, the ManaMed logo, are trademarks of ManaMed LLC
Registered in the U.S.