

CLOSED**EN**

UNDERLAY TT

SIZING GUIDE

Smart Size Guide

[Click here](#)


Intended for all suspensions systems:

- Distal locking (pin or lanyard)
- Active / passive vacuum systems
- Suspension sleeve

To determine the most suitable size for the patient, measure the circumference at 10 cm from distal end. Then measure the length of the trimmed liner.


STEP 1

Which size?

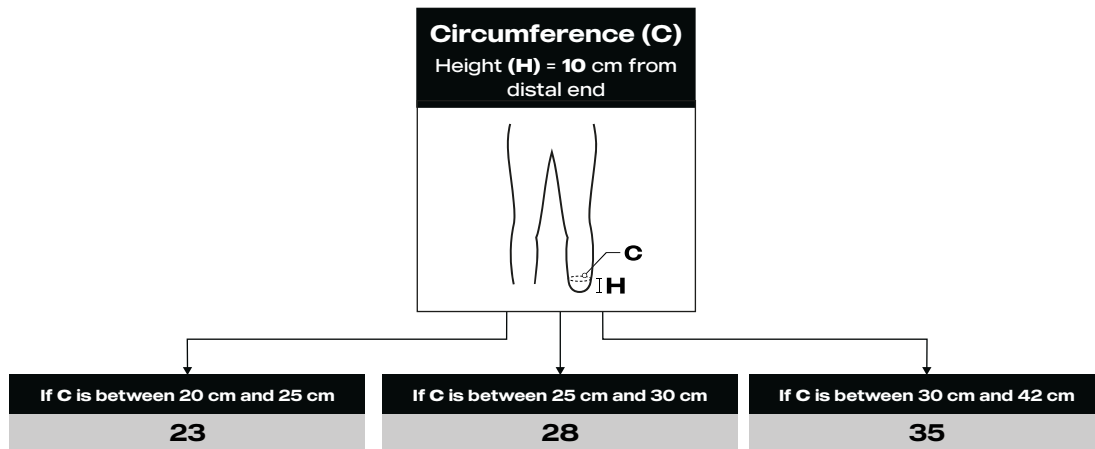

STEP 2

Which length?

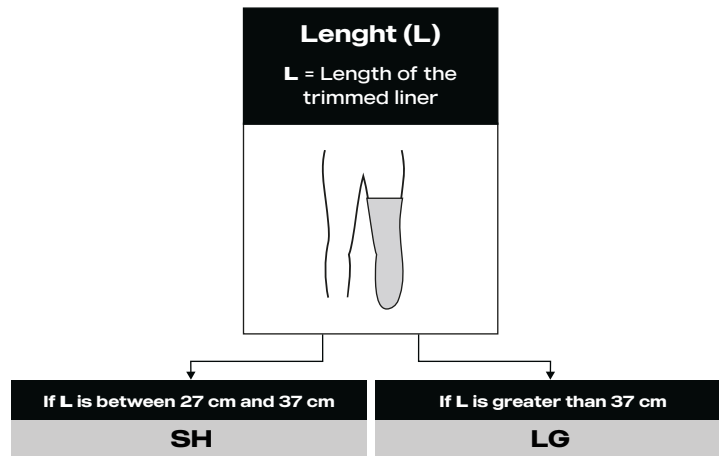


Circumference (C)	Length (L)	Code
20-25 cm	27-37 cm	UDTT-23-SH-C
	37+ cm	UDTT-23-LG-C
25-30 cm	27-37 cm	UDTT-28-SH-C
	37+ cm	UDTT-28-LG-C
30-42 cm	27-37 cm	UDTT-35-SH-C
	37+ cm	UDTT-35-LG-C


Intended for:




Pin

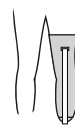


Lanyard

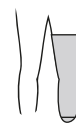


Cushion

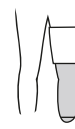

Suspension sleeve

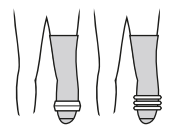

Seal-in liner

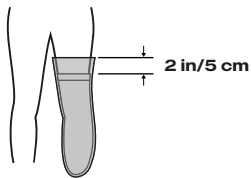Important:Ensure the Underlay Closed sits **2 in/5 cm** under the top edge of the liner.

See next page for
Underlay TT Open →

ETHNOCARE™



[Click here](#)

Intended for distal locking systems:
 -Pin system
 -Lanyard system
 -Distal seal

To determine the most suitable size for the patient, measure two circumferences over the skin: at 10 cm from distal end and just above the trimmed liner. Then measure the length of the trimmed liner.


STEP 1

Which size?

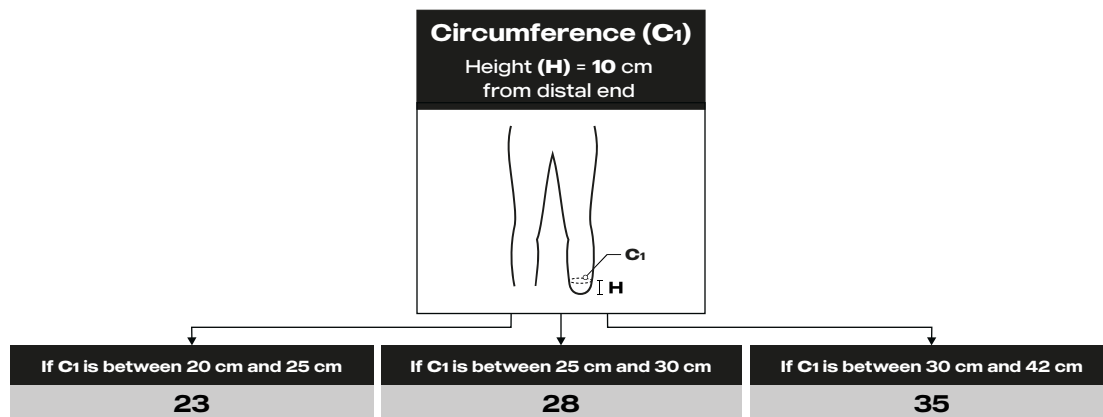


STEP 2

Which version?

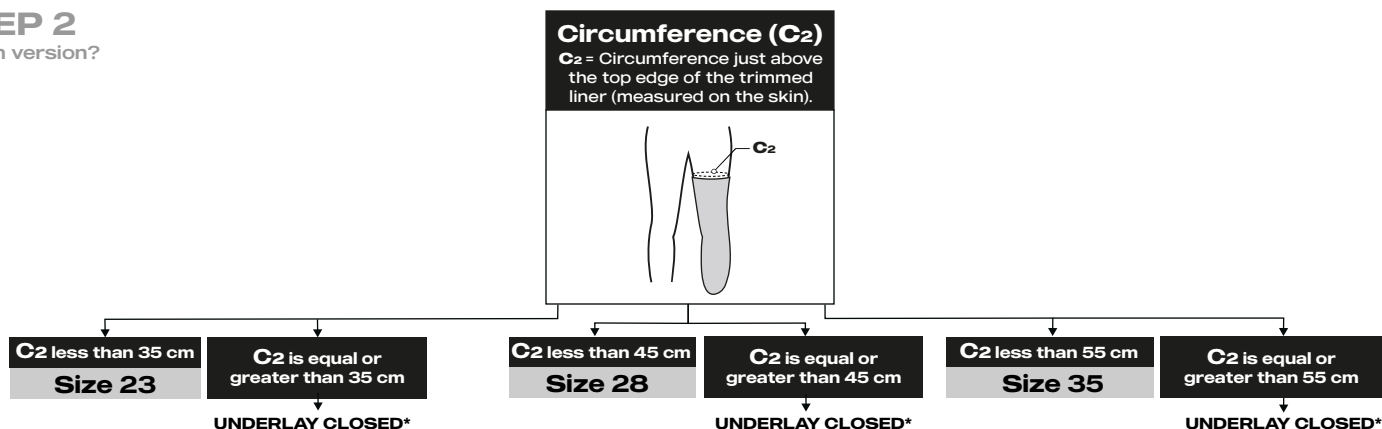


STEP 3

Which length?

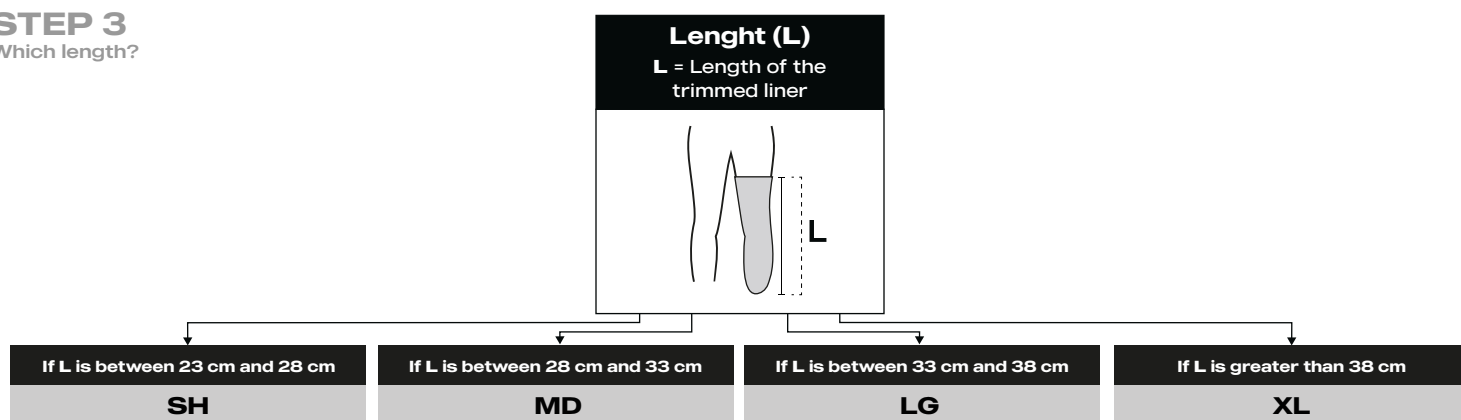

Circumference (C ₁)	Length (L)	Code
20-25 cm	23-28 cm	UDTT-23-SH
	28-33 cm	UDTT-23-MD
	33-38 cm	UDTT-23-LG
	38+ cm	UDTT-23-XL
25-30 cm	23-28 cm	UDTT-28-SH
	28-33 cm	UDTT-28-MD
	33-38 cm	UDTT-28-LG
	38+ cm	UDTT-28-XL
30-42 cm	23-28 cm	UDTT-35-SH
	28-33 cm	UDTT-35-MD
	33-38 cm	UDTT-35-LG
	38+ cm	UDTT-35-XL


Intended for:

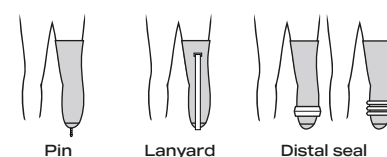

*If the measurement exceeds the indicated C₂ value, the limb is too conical for available sizes. Select an **Underlay Closed** for optimal fit, as it sits completely under the liner.

Important:

Ensure the Underlay Open sits 4 in/10 cm from distal end.

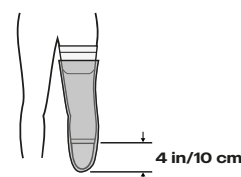

See next page for
Underlay TT Closed →

ETHNOCARE™